

Create it 

TIP GUIDE

AIR SCALER

NSK

A Tip range to suit
your needs.


Scaling P04



Perio P06



Perio (Maintenance Tip) P07-08

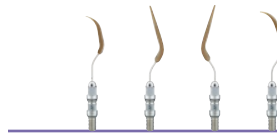


Restorative (For Minimal Intervention) P10




Restorative (For Preparation of crown cores) P12-14




V-TIP SYSTEM (Endodontics) P16

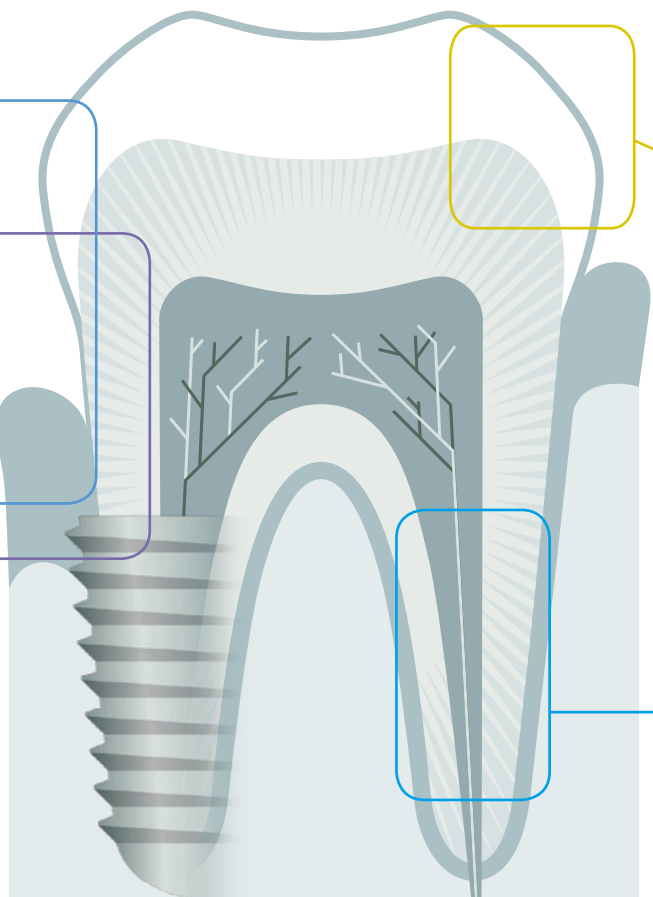

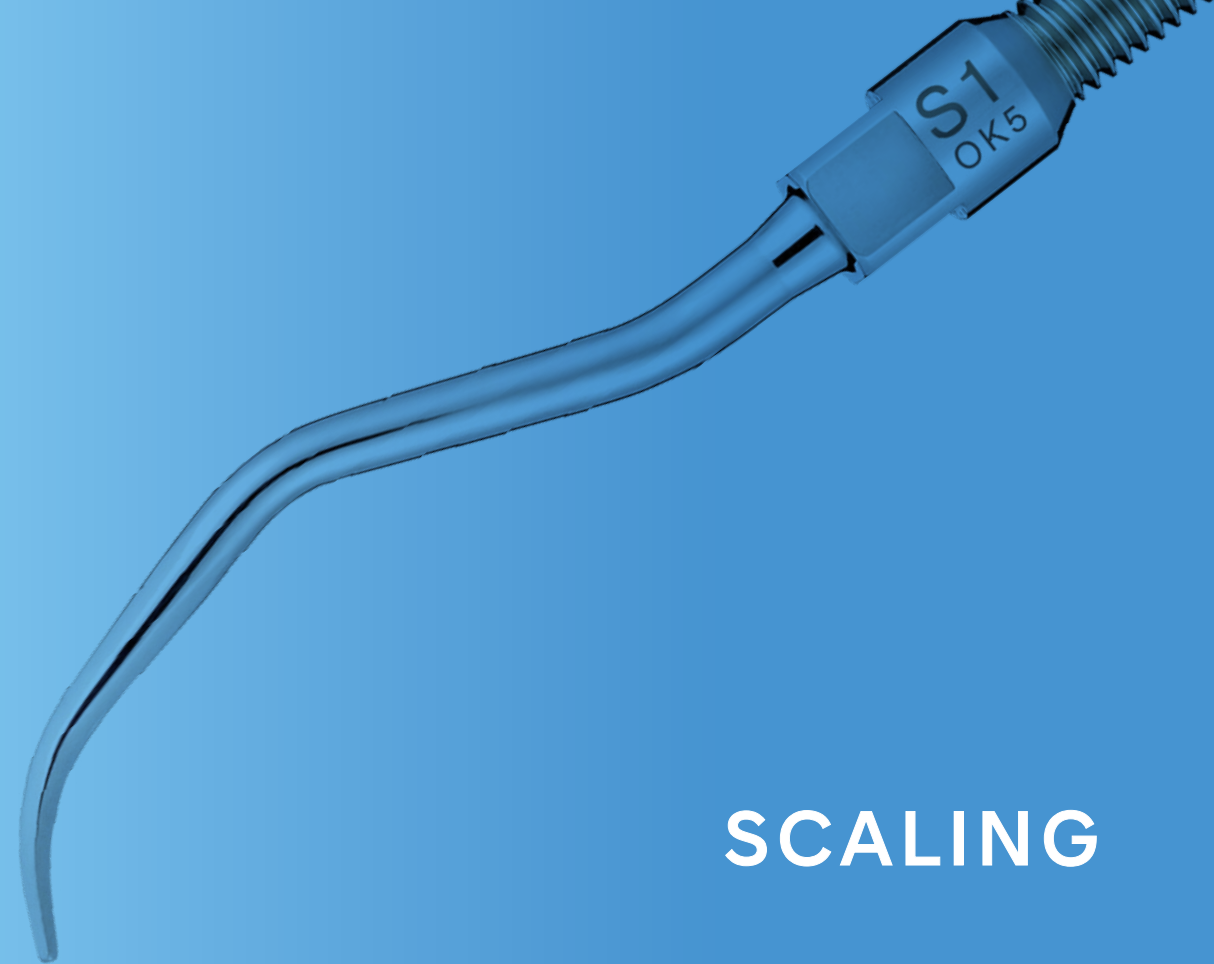



SCALING

SCALING

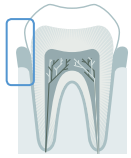

MODEL **S1**

- Removal of Supra gingival calculus.
- Removal of marginal gingival calculus.

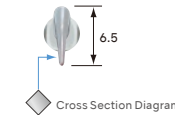
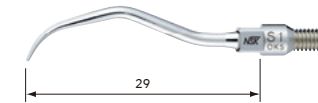


ORDER CODE **Z252411**

- Universal Type

MODEL **S2**

- Removal of Supra gingival calculus.

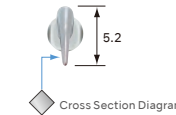
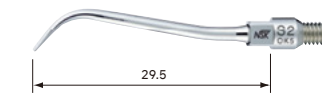


ORDER CODE **Z252413**

- Sickle Type


MODEL **S3**

- Removal of Supra gingival calculus.
- Removal of marginal gingival calculus.

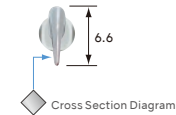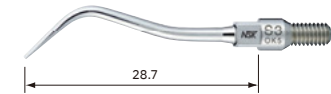

ORDER CODE **Z252412**


- Perio Type

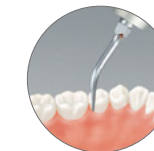





PERIO

PERIO

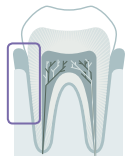

MODEL S10

- Removal of sub gingival calculus.
- Removal of narrow interdental space calculus.

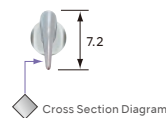
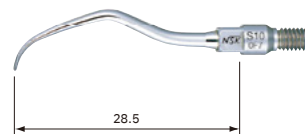


ORDER CODE **Z252420**

- Universal Type Slim Line

MODEL S20

- Removal of sub gingival calculus.
- Removal of sub gingival biofilm

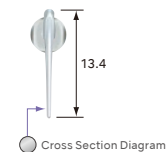
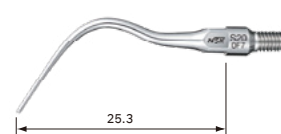


ORDER CODE **Z252441**

- Root Planing (Straight tip)

MODEL S21R

- Removal of sub gingival calculus.
- Removal of narrow interdental space calculus.
- Removal of sub gingival biofilm

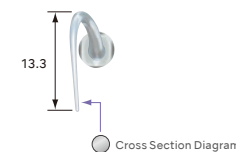
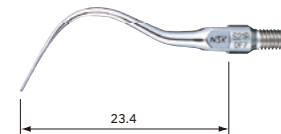


ORDER CODE **Z252442**

- Root Planing (Right curved tip)

MODEL S21L

- Removal of sub gingival biofilm

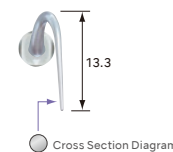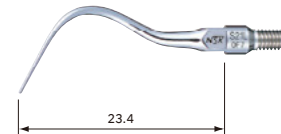


ORDER CODE **Z252443**

- Root Planing (Left curved type)




PERIO (Maintenance Tip)

For the maintenance of prostheses, implant-supported prostheses.
These tips are made from a special material which will not damage prostheses and restorative materials.

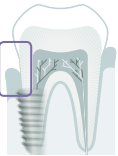

MODEL **V-P10**

V-P10 is for the maintenance of prostheses, its shape offers easy access to all treatment areas.

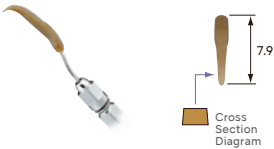


ORDER CODE **Y900184**

• Pack of 3 • S35 holder is not included



MODEL **V-P11R**

V-P11R/L left and right are double bent tips. This tip is extra long to easily access difficult to reach areas and to avoid interference with neighbouring teeth or contours. Its shape was designed to best fit different tooth morphologies, making this a versatile tip for a wide range of uses from anterior teeth to premolars.

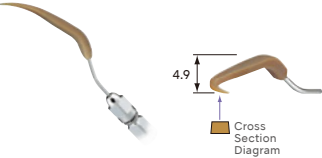

ORDER CODE **Y1002165**

• Right curved type • Pack of 3 • S35 holder is not included



MODEL **V-P11L**




ORDER CODE **Y1002166**

• Left curved type • Pack of 3 • S35 holder is not included



MODEL **V-P12**


V-P12 is a longer version of the popular V-P10 tip. It serves a wide range of applications.

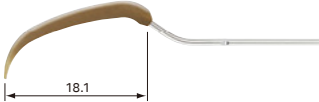

ORDER CODE **Y1002167**

• Pack of 3 • S35 holder is not included


MODEL **S35**

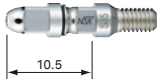

ORDER CODE **Z252015**

• Includes E-Tip replacement wrench • Plastic Tip is not included


MODEL **Perio-Control Tips Air Scaler Kit**

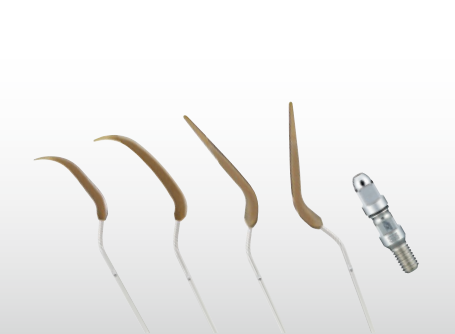


ORDER CODE **Y1002754**

• S35 holder : 1pc.
• Includes E-Tip replacement wrench
• V-P10, V-P11R, V-P11L, V-P12 : 1 pc. each

RESTORATIVE





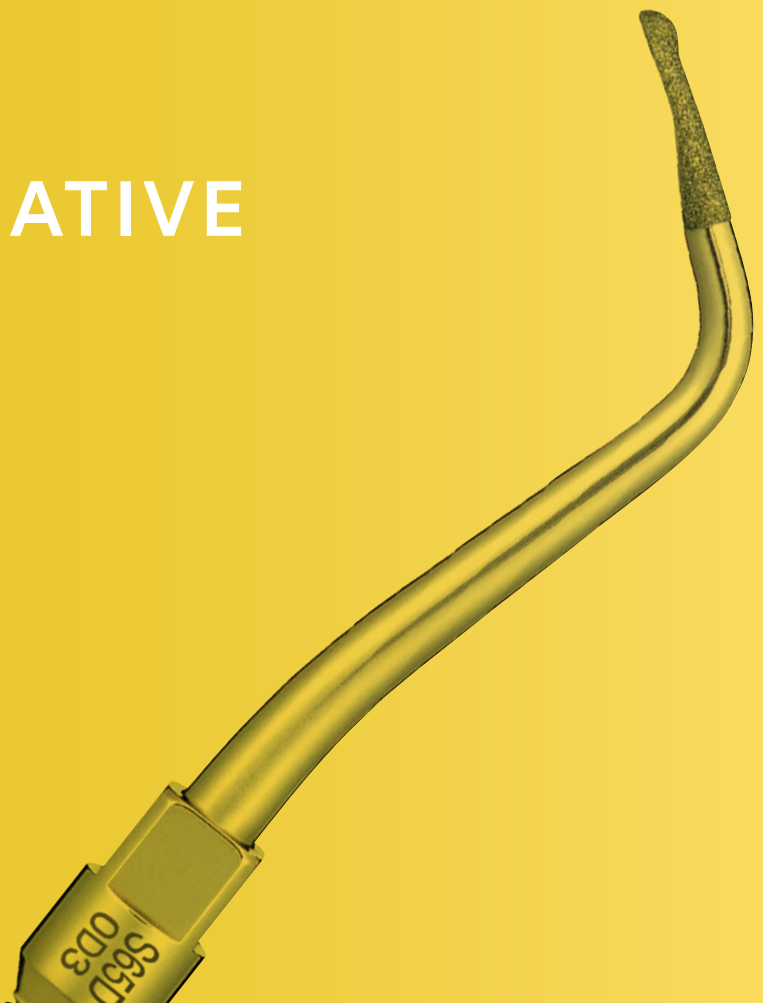


RESTORATIVE (For Minimal Intervention)

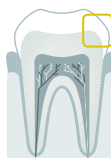

| | | | | | | | |
|-------|-------------|-------|-------------|-------|-------------|-------|-------------|
| MODEL | S65D | MODEL | S66D | MODEL | S67D | MODEL | S68D |
|-------|-------------|-------|-------------|-------|-------------|-------|-------------|

The spherical shape of these tips make them suitable for delicate adjustments as well as caries removal in the cervical and interproximal region.

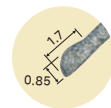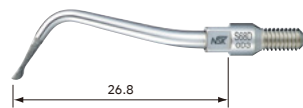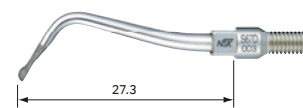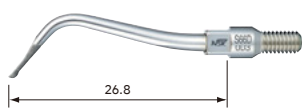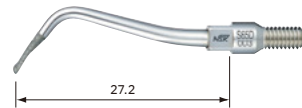

ORDER CODE **Z252465**

• Diamond coated

ORDER CODE **Z252466**

• Diamond coated

ORDER CODE **Z252467**

• Diamond coated


ORDER CODE **Z252468**

• Diamond coated



Caries removal in the cervical and interproximal region.

RESTORATIVE

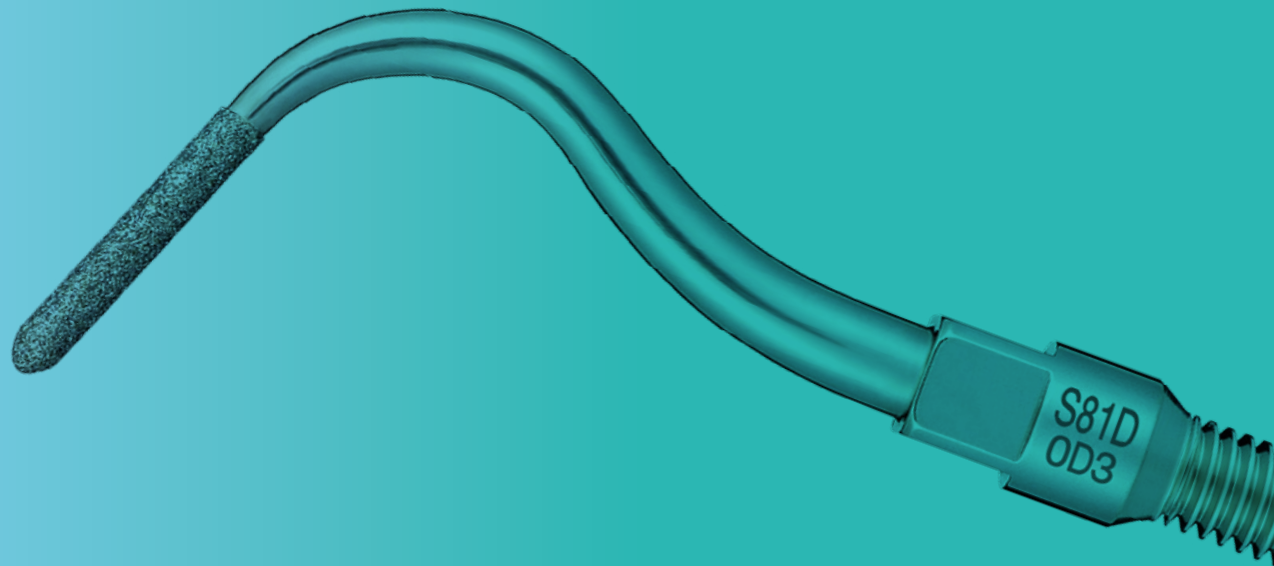



RESTORATIVE (For Preparation of crown cores)



MODEL **S71D**

MODEL **S72D**

These tips are safe sided and flat for minimal mesial or distal preparations.

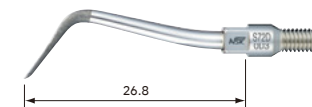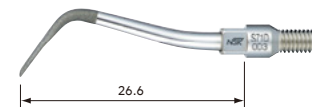


ORDER CODE **Z252471**

ORDER CODE **Z252472**

• Diamond coated

• Diamond coated


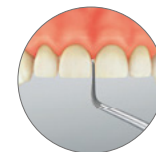


RESTORATIVE (For Preparation of crown cores)




MODEL **S81D** MODEL **S82D** MODEL **S83D**

The tips can be used for finishing margins after prepping the area with a conventional dental bur.
This will enable you to finish minute sharp edges and gingival sulcus without damage, something that cannot be accomplished with rotary instruments alone.

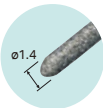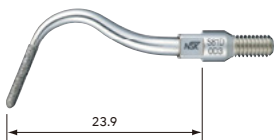




ORDER CODE **Z252469**

• Diamond coated

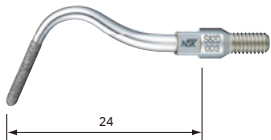



ORDER CODE **Z252462**

• Diamond coated

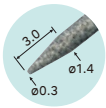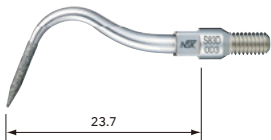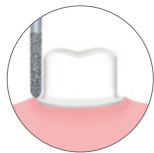


ORDER CODE **Z252461**


• Diamond coated

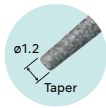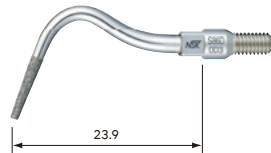

ORDER CODE **Z252460**

• Diamond coated

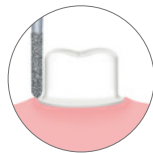


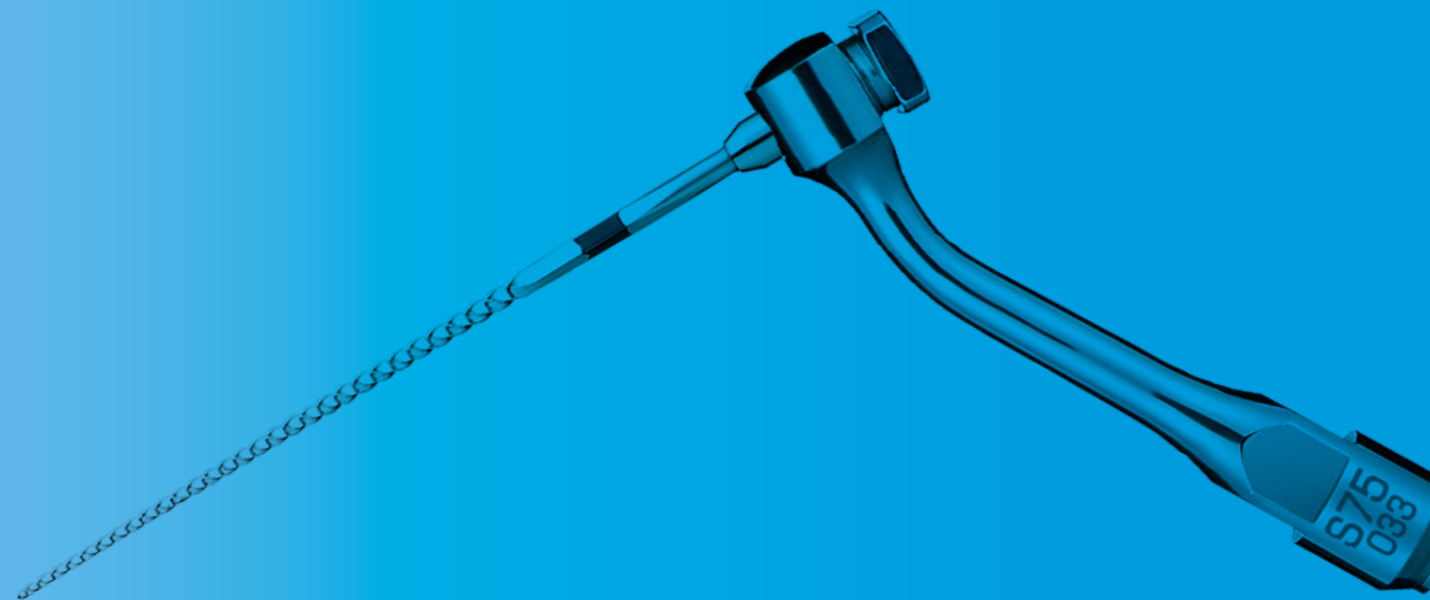

V-TIP SYSTEM

V-TIP SYSTEM

NSK V-Tip system is innovative, multi-functional and economical. The V-Tip holder is a separate component to the actual cutting tips and can be easily and inexpensively replaced. NSK V-Tip holder system offers the flexibility and expandability to meet the requirements for multi-purpose clinical applications.

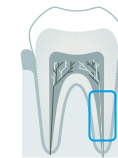

MODEL **S75**



ORDER CODE **Z252035**


• Included E-Tip replacement wrench

V-TIP SYSTEM (Endodontics)



For Root Canal Cleaning


MODEL **V-U15**

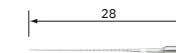

ORDER CODE **Y900177**

• V-file Irrigation 23 mm • ISO 15 • Pack of 6


MODEL **V-U20**

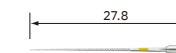

ORDER CODE **Y900178**

• V-file Irrigation 23 mm • ISO 20 • Pack of 6


MODEL **V-U25**

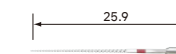

ORDER CODE **Y900179**

• V-file Irrigation 23 mm • ISO 25 • Pack of 6


MODEL **V-U30**

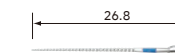

ORDER CODE **Y900180**

• V-file Irrigation 23 mm • ISO 30 • Pack of 6

MODEL **V-U35**

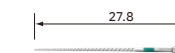

ORDER CODE **Y900181**

• V-file Irrigation 23 mm • ISO 35 • Pack of 6



Ti-Max S970 Air Scaler



- Titanium Body with Scratch Resistant DURACOAT
- 3 level Power Ring • Cellular Glass Optics
- Includes 3 tips (S1,S2 & S3), Tip Cover(L) and Tip Wrench(CR-20)

NSK Couplings

| Optic | MODEL | ORDER CODE |
|-------|-------|------------|
| ● | S970L | T1005 |
| — | S970 | T1020 |

KaVo® Coupling

| Optic | MODEL | ORDER CODE |
|-------|--------|------------|
| ● | S970KL | T1013 |

Dentsply Sirona® Coupling

| Optic | MODEL | ORDER CODE |
|-------|--------|------------|
| ● | S970SL | T1015 |

KaVo and Multiflex are registered trademarks of KaVo Dental GmbH.
Dentsply Sirona and Sirona are registered trademarks of Dentsply Sirona Inc.

Washable in the thermodisinfectior Autoclavable

POWER GUIDE

NSK Power Ring offers three different power level settings for diverse clinical needs.
While the gentle minimum mode is for delicate procedures the maximum power mode is good for general scaling.

| APPLICATIONS | TIPS | POWER RING SETTING | | |
|--|---|--------------------|---|---|
| | | 1 | 2 | 3 |
| Scaling | S1, S2, S3 | ←→ | | |
| Perio | S10, S20, S21R, S21L | ←→ | | |
| Perio (Maintenance Tip) | S35+ V-Tip (V-P10, V-P11R, V-P11L, V-P12) | ←→ | | |
| Restorative (For Preparation of crown cores) | S81D, S82D, S83D, S86D, S71D, S72D | ←→ | | |
| Restorative (For Minimal Intervention) | S65D, S66D, S67D, S68D | ←→ | | |
| V-TIP SYSTEM (Endodontics) | S75 + V-Tip (V-U15~V-U35) | ←→ | | |

- Always start out with lower power setting and then adjust to higher setting if necessary.
- Incorrect operation may cause tip breakage even if tip is used with in the recommended power range. Consult the operating manual before use.

NAKANISHI INC. www.nsk-dental.com
700 Shimohinata, Kanuma, Tochigi 322-8666, Japan

NSK Europe GmbH
Elly-Beinhorn-Str.8
65760 Eschborn
Germany

NSK United Kingdom Ltd.
Hertford House
Rutherford Close, Stevenage
Hertfordshire, SG1 2EF, UK

NSK Dental Nordic
Furstenbergsgatan 4
416 64 Gothenburg
Sweden

NSK Oceania Pty. Ltd.
Unit 12, 809-821 Botany Road
Rosebery NSW 2018
Australia

NSK Asia
No. 8, T-ONE Building, Room No. 26-109
26th Floor, Soi Sukhumvit 40
Sukhumvit Road, Bangkok 10110 Thailand

NSK Middle East FZCO
Room 6EA701, 7th Floor
East Wing No.6
Dubai Airport Free Zone
PO Box 54316, Dubai, UAE

Specifications are subject to change without notice.