





NSK Dental LLC

www.nskdental.com

1800 Global Parkway, Hoffman Estates, IL 60192-1578, USA





Vario Surg Handpieces Optic / Non-Optic

A Powerful Ultrasonic Surgical System with LED

Applications

- Bone surgery
- Sinus lifts
- Periodontal surgery
- Endodontic surgery

LED

Twin LED lights



NSK LED delivers natural daylight quality light bringing significant advantages to the way you work

Clearer: LEDs generate natural daylight quality light which assures clearer vision when compared to halogen light.

Durable: LED light sources are longer lasting and more durable than halogen bulbs.

Comfortable Grip Handpiece

Superior ergonomics of the VarioSurg handpiece, including the slim and lightweight body, ensures minimal hand fatigue and greater accessibility to difficult oral regions. The twin LEDs integrated at the head allows for outstanding visibility, thus treatment is safer and precisely accurate. Oscillation at the tip of the handpiece due to advanced ultrasonics eliminates heat at the operation site, allowing for minimally invasive procedures.

Shorter and Lighter Handpiece

A shorter and lighter NSK LED handpiece is achieves a perfect balance. Surgical procedures are easier and stress free to help you provide more precise and efficient treatment.



The Special Tip of VarioSurg made by NSK's High Construction Accuracy

Maximizing sharpness of tips through the three-dimensionally cutting (toothing) to the blade edge.



3 VarioSurg

The Advanced Ultrasonic Technology of NSK

VarioSurg is equipped with various advanced functions of ultrasonic bone surgery through it's specially-made engine. The Feedback & Auto Tuning Function achieve high cutting efficiency through the automatic power control to exacting levels according to grinding surface.

The Burst Function generates the "hammer effect", which adds more powerful vibration to the vibrating tip, thus improving the ability to cut hard tissue.

Advantages

Power: superior cutting capability shortens treatment time

Minimally invasive: eliminates heat generation during bone dissection

Safe: avoids damage to soft tissue even during accidental contact

Precision: cavitations effect of irrigation solution keeps surgical site blood free

Burst Mode

Burst Mode creats a hammer drill effect promoting cutting of the hardest tissues. Superior cutting capability of the burst mode facilitates cutting of hard tissue without stress, and effectively shortens treatment time.

Continuous vibration mode

Burst mode

Feedback Function

Feedback Function constantly checks the performance parameters of the ultrasonic unit during operation.


Simultaneously it controls the power output, which enables to supply the power suited for bone cutting procedure.

Simple, intuitive graphic operation panel



Main Features

- 3 modes (SURG, ENDO, PERIO)
- Burst Mode
- 9 programmes (SURG x 5, ENDO x 2, PERIO x 2)
- Automatic Cleaning Mode
- Ergonomic Handpiece (LED or Non-Optic)
- Feedback Function

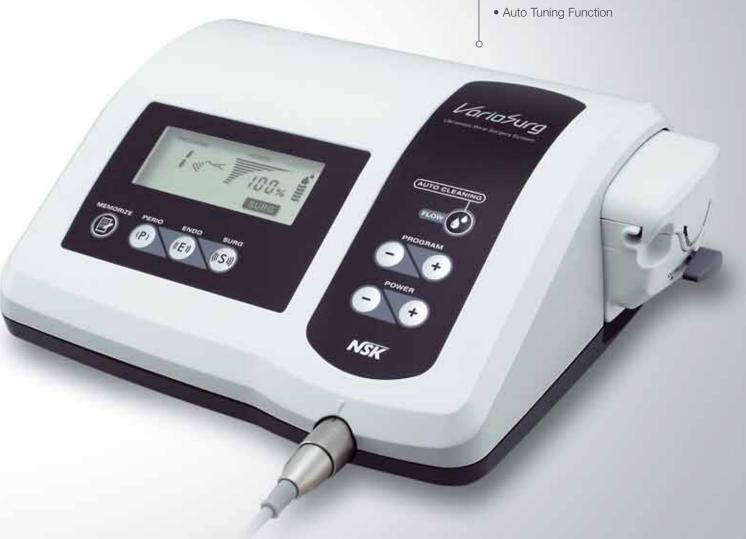

5 VarioSurg VarioSurg

Variosurg Complete Sets Optic / Non-Optic

VarioSurg Complete Set with LED Handpiece

MODEL: VSRG OPT ORDER CODE: Y1001928

Contents

- Control Unit LED Handpiece with 2m cord
- FC-51 Foot Control Sterilisation Cassette
- Handpiece Stand Irrigation tube (5 pcs.) and other accessories
- Basic-S Kit (SG1, SG3, SG5, SG6D, SG7D, SG11 and Tip holder)

VarioSurg Complete Set with Non-Optic Handpiece

MODEL: **VSRG 120V2** ORDER CODE: **Y1001939**

Contents

- Control Unit Non-Optic Handpiece with 2m cord
- FC-51 Foot Control Sterilisation Cassette
- Handpiece Stand Irrigation tube (5 pcs.) and other accessories
- Basic-S Kit (SG1, SG3, SG5, SG6D, SG7D, SG11 and Tip holder)



Specifications

Control Unit

- Frequency 28-32 kHz
- Power Supply : AC120V 50/60 Hz
- Irrigation Flow Rate : 9~90 mL / min
- Programs SURG x 5, ENDO x 2, PERIO x 2 : W268 x D230 x H103 (mm)
- Dimensions

Foot Control


- Foot Control Functions :
- ON-OFF Switch for Water Flow.
- PROGRAM (+) Switch, PROGRAM (-) Switch

Vibration ON-OFF Switch

Foot Control

Hands-free Program Adjustments Via Foot Control

The VarioSurg Foot Control, which meets IPX8 Standards for medical foot control systems, allows perfect control during delicate surgical procedures. A bracket enables the Foot Control to be easily raised or moved without using your hands during use.

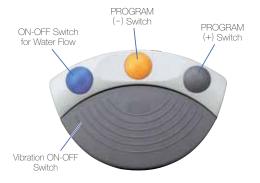



MODEL: FC-51 ORDER CODE: Z316001

• With 2m cord

Sterilization Cassette

The Sterilization Cassette is designed for the safe processing and storage of VarioSurg components.





MODEL: VA-SG-CASE ORDER CODE: 20001326

- Dimensions : W281 x D171.5 x H47 (mm)
- With appropriate compartments for handpiece, cord, tip replacement wrench and tip holders

Carrying Case (Optional)

The NSK Carrying Case can accommodate all VarioSurg components.



MODEL: Carrying Case ORDER CODE: Y1001907 (VarioSurg)

• Dimensions : W469 x D344 x H180 (mm)

7 VarioSurg

VarioSurg Ultrasonic Tips



MODEL: SG14L ORDER CODE: Z305123 • TiN coating • Left curved tip*

• Five teeth; 0.6mm thick



Unilateral anterior mandibular sampling



Excision of bone torus



Anterior mandibylar harvesting



Mandibular harvesting at approach stage





MODEL: SG3 ORDER CODE: **Z305104**

• TiN coating • Triple sided edge spatula tip

• TiN coating • Sharp edge • 0.5mm thick



MODEL: SG4 ORDER CODE: Z305105

• TiN coating • Edge spatula tip



Ostroplasties



MODEL: SG5 ORDER CODE: Z305106

• TiN coating • Rounded edge spatula tip



Avulsion of wisdom teeth and impacted teeth



Sinus Lift



MODEL: SG6D ORDER CODE: Z305107

Diamond coating



MODEL: SG7D ORDER CODE: Z305108

· Diamond coating

Sinus graft





MODEL: **SG17** ORDER CODE: **Z305132**

• TiN coating • 0.7mm thick

0.7mm thick



MODEL: SG18R ORDER CODE: Z305133 • TiN coating • Right angled tip* 0.7mm thick

MODEL: SG18L ORDER CODE: Z305134 • TiN coating • Left angled tip*

Extraction

Sinus Membrane Tool



MODEL: **SG9** ORDER CODE: **Z305110**

• Flat circular convex elevator • Angled at 90°



MODEL: **SG10** ORDER CODE: **Z305111**

• Flat circular convex elevator • Angled at 135°




MODEL: **SG11** ORDER CODE: **Z305112**

• Cone compressor



Sinus membrane

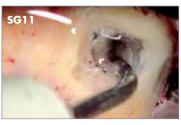

Sinus graft









MODEL: SG16A ORDER CODE: Z305128

• Diamond coating • Diameter; 2.0mm



MODEL: SG16B ORDER CODE: Z305129

• Diamond coating • Diameter; 2.6mm



Conical preparations for implant osteotomies




Implant preparation

Scaling



Perio (Root Planing)

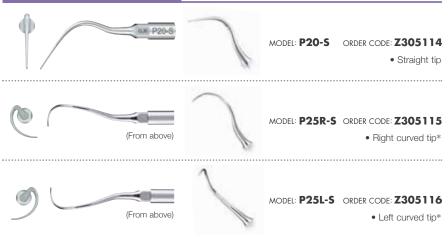

Maintenance (V-Tip)



Retrograde Endo



MODEL: **E32D-S** ORDER CODE: **Z305121**

• Diamond coated

• For anterior teeth (90°)

Tip Kits



MODEL: Basic-S Kit ORDER CODE: **Y900687**

- SG1, SG3, SG5, SG6D, SG7D, SG11
- Tip holder



MODEL: Bone Cut Kit ORDER CODE: **Y900688**

Contents

- SG1, SG2R, SG4, SG2L, SG6D
- Tip holder



MODEL: Sinus Lift Kit ORDER CODE: **Y900689**

Contents

- SG1, SG3, SG6D, SG9, SG10, SG11
- Tip holder



MODEL: Endo-S Kit ORDER CODE: **Y900691**

Contents

- G1-S, E30RD-S, E30LD-S, E31D-S, E32D-S
- Tip holder



MODEL: Implant Preparation Kit ORDER CODE: **Y900774**

Contents

- SG15A, SG15B, SG16A, SG16B
- Tip holder

Spare Parts

	MODEL	ORDER CODE
LED X E	VS-LED-HPSC LED Handpiece (with Cord 2m)	E1084
一	VA-SGN-HPSC Non-Optic Handpiece (with Cord 2m)	E377
	E Tip Replacement Wrench For V10-S	Z217399
	VS Tip Wrench For SG1A, SG8A (CR-40)	Z305350
X EE	Tip Wrench (CR-30)	10000977

	MODEL	ORDER CODE
X W	Tips Holder	20001327
A W	Handpiece Stand	Z1029201
	Irrigation Tube (Pack of 5)	Y900113
X III	Tube Holder (Pack of 7)	Y900767
X III	VA-SG-CASE Sterilization Cassette	20001326